

Designed to hole cubed optics or right angle prismms.

All four polished surfaces of a cubed prism can be used with these mounts.

- The KKD series has tip, tilt and yaw rotation for fine adjustment of reflected beams.
- The clamp screw will not rotate the prism when clamping.



Guide

- ► θαβ axis stages without prism retainer and strut (KKD-C) are also available. ► WEE Reference Catalog Code W7111
- ▶ Post length can be changed by specifying the post length when you place an order. We may charge the difference in price depending on the length. Contact our Sales Division for more information.


Attention


- ▶ KKD-25PH and PLH-25 might not exert sufficient pressure to retain a prism, and there is a risk that optics might fall out. Use after making sure that the prism is fixed.
- ▶ After adjusting the KKD series, if the prism retainer is lifted, the prism table sometimes moves, throwing off adjustment of tilt and rotation. Use without touching the prism retainer after adjustment.




Outline Drawing KKD-25PH Pan head screw M2×6...3 screws Rotational adjustment 35 3 20 Tilt adjustment 12.5 Prism retainer



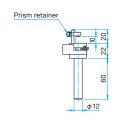

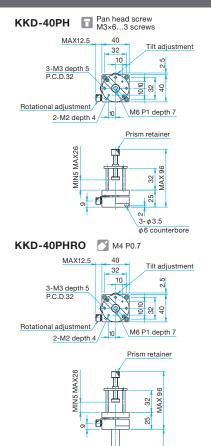

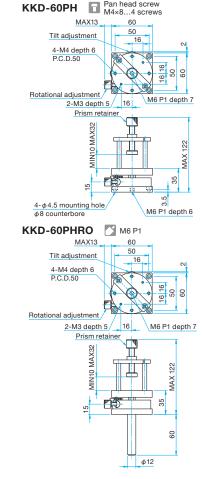

With Tilt and Rotational Adjustment						
Part Number	Compatible Optics Dimensions [mm]	Adjustment Range Rotation [°]	Adjustment Range Tilt [°]	Reso Rotation [°/rotation]	lution Tilt [°/rotation]	Weight [kg]
KKD-25PH	- □10	±3	±3	about 0.9	about 2.2	0.05
KKD-25PHRO	– □10	±3	±3	about 0.9	about 2.2	0.11
KKD-40PH	□5 – □26	±3	±3	about 2.0	about 1.5	0.20
KKD-40PHRO	□5 – □26	±3	±3	about 2.0	about 1.5	0.25
KKD-60PH	□10 – □32	±3	±3	about 1.7	about 1.0	0.40
KKD-60PHRO	□10 – □32	±3	±3	about 1.7	about 1.0	0.45

φ12

Application Systems

Optics & Optical Coatings

Opto-Mechanics

Bases

Manual Stages

Actuators & Adjusters

Motoeized Stages

Light Sources & Laser Safety

Index

Guide

Mirrors Lenses

Prisms

Polarizers

Lasers

Beam Shaping Diffusers

Filters

Shutter

Others

Fiber



Prism Holders

KKD/PLH/PAD



Application Systems

Optics & Optical Coatings

Opto-Mechanics

Bases

Manual **Stages**

Actuators & Adjusters

Motoeized **Stages**

Light Sources & Laser Safety

Index

Guide

Mirrors Lenses

Prisms

Polarizers

Lasers

Beam Shaping Diffusers

Filters

Shutter

Others Fiber

Outline Drawing


Prism retainer

PLH-25 M4 P0.7








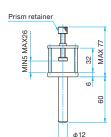

PLH-60 M6 P1



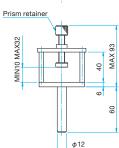

Specifications	Primary material: Aluminum Finish: Black Anodized		
Part Number	Options specified*	Compatible Optics Dimensions [mm]	Weight [kg]
PLH-25	N	- □10	0.08
PLH-40	N	□5 – □26	0.14
PLH-60	N	□10 – □32	0.29

^{*} For specifying options, please refer to "Conversion of Posts, Post Holders and Pedestal Bases of Holders". Reference) C007

Prism Holder








Compact Mount for cube beamsplitters. By clamping on the lower edges of the cube more of the prism faces are available for use.

• It is possible to put close to other optical elements because the holder shapes match to the width of the prism.

Outline Drawing

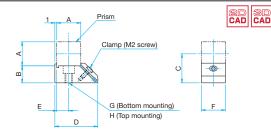

Part Number	A (mm)	B (mm)	C (mm)	D (mm)	E (mm)	F (mm)	G	Н
PAD-10	10	10	15	20	5	10	M4	МЗ
PAD-12.7	12.7	8.65	15	22.7	6.35	12.7	M4	МЗ
PAD-15	15	12.5	20	25	7.5	15	M4	МЗ
PAD-20	20	15	25	35	10	20	M6	M4
PAD-25.4	25.4	12.3	25	39	12.5	25.4	M6	M4

Guide

▶ Can be mounted on posts, pillars or directly on a baseplate or stage with an M4 thread.

Attention

- ▶ Over tightening the clamp may break glass. Please tighten with the minimum necessary force to insure the cube does not move.
- Recommended for use with prisms with outer dimension tolerance of ±0.2mm.

		naterial: Aluminum ack Anodized
Part Number	Compatible Optics Dimensions [mm]	Weight [kg]
PAD-10	□10	0.005
PAD-12.7	□12.7	0.006
PAD-15	□15	0.01
PAD-20	□20	0.022
PAD-25.4	□25.4	0.026