

Vacuum Mirror Holder | VMHG

NEW

RoHS

Catalog Code

W4150

Low out-gas mirror holders for installing optical system in the high vacuum chamber.

- These can be used for UV laser optical system and in the clean room.
- **Adjustment screw is removable (except for VMHG-12.7). This enables installing electric actuators of remote control device. (Actuators are not included.)**
- Screw hole grooving and vent hole processing are applied to prevent the out-gas emits when forming vacuum.
- BARRIERTA SUPER IS/V is used for the grease of sliding part.



Guide


- ▶ We deliver the item with specialized vacuum packaging to keep its high clean state from the factory to the destination.
- ▶ Please use vacuum screws (VPCB/VAGB) for installing base plates. [WEB Reference](#) **Catalog Code** W6111
- ▶ We provide vacuum pedestal bases (VPST) and posts (VRO). [WEB Reference](#) **Catalog Code** W6108, [WEB Reference](#) **Catalog Code** W6107
- ▶ We will produce custom products for the requested specifications.

Attention


- ▶ Specialized nut turning tool (VNS-12) is necessary in removing adjustment screw or equipping actuators. W6112
- ▶ Please open the package in the clean environment. Please do not touch it by your bare hands.

Outline Drawing


VMHG-12.7

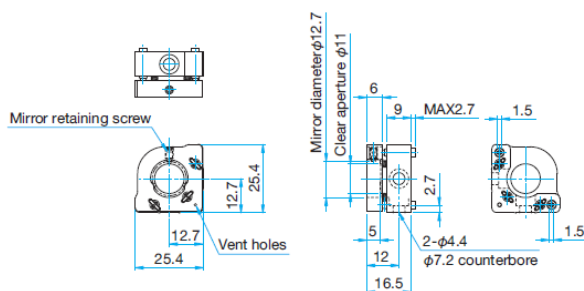


VMHG-25.4

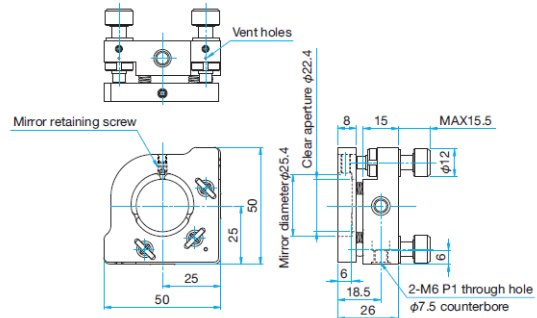

VMHG-30

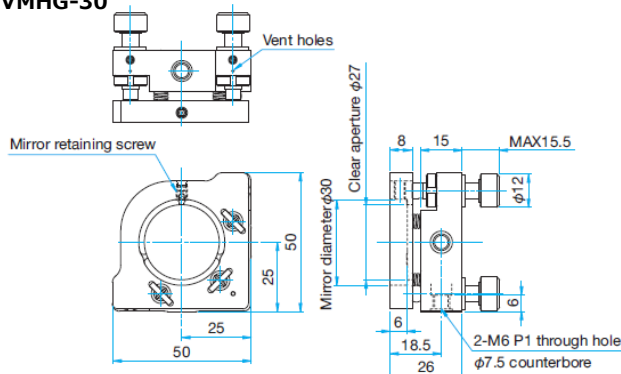

VMHG-50.8

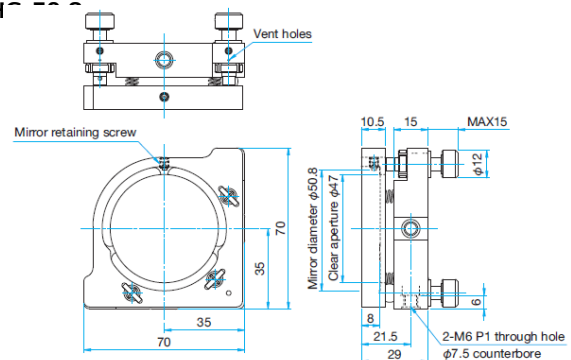

Specifications										Primary material: Steel (SUS316L) Finish: None
Part Number	Compatible Optics Diameter [mm]	Compatible Optics Thickness [mm]	Clear Aperture [mm]	Number of Adjustment Axes	Adjustment Range		Resolution		Weight [kg]	Compatible Mounting Screw Number
VMHG-12.7	φ12.7	3~7	φ11	2	Tilt [°]	Rotation [°]	Tilt [°/Rotation]	Rotation [°/Rotation]	0.032	VPCB-4-8 or VAGB-4-8
VMHG-25.4	φ25.4, φ25	3~10	φ22.4	2	±3	±3	about0.39	about0.39	0.21	VPCB-4-10 or VAGB-4-10
VMHG-30	φ30	3~10	φ27	2	±3	±3	about0.39	about0.39	0.19	VPCB-4-10 or VAGB-4-10
VMHG-50.8	φ50, φ50.8	3~10	φ47	2	±2.5	±2.5	about0.26	about0.26	0.35	VPCB-4-10 or VAGB-4-10



SIGMAKOKI

Tokyo Head Office
5F, SIGMA KOKI Tokyo Head office 1-19-9, Midori, Sumida-ku, Tokyo, 130-0021, JAPAN
TEL. +81-3-5638-8228
FAX. +81-3-5638-6550
E-mail international@sigma-koki.com

Osaka Branch
4-9-28 Nishi-Nakajima, Yodogawa-ku, Osaka 532-0011
TEL. +81-6-6307-4835
FAX. +81-6-6307-4834
E-mail sales.osaka@sigma-koki.com

Kyushu Sales Office
3-17 Hie-machi, Hakata-ku, Fukuoka 812-0014
TEL. +81-92-481-4300
FAX. +81-92-481-4310
E-mail sales.kyushu@sigma-koki.com